

AVTECH

High Performance FREEZE DRYING SYSTEMS



LyoPro series Laboratory freeze dryer

Optimize and simplify your freeze-drying process

Antech Scientific offers advanced LyoPro series freeze-drying system with unique LyoSmart operating system, realizing full control of vacuum degree, making the freeze-drying process simpler, faster and better.

LyoPro series freeze dryers are compact in design, high-efficiency in performance, smart in operation, and in daily use.



► Innovative features

- Compact, high-performance laboratory systems with a small footprint
- 7" touch screen installed with LyoSmart operating system
- Programmable freeze-drying steps
- Controlled vacuum degree for optimized process times
- Built-in ice condenser coils for higher efficiency
- Fast defrost via integrated hot gas
- Automatic non-return venting valve, prevents "vacuum pump" oil mist from entering cold trap
- Pirani vacuum sensor with accuracy 0.001 mbar
- Low noise level: 48~52dB(A)
- High quality vacuum pump
- PTFE coated stainless steel ice condenser

► LyoSmart operating system

7" touch screen controller, installed with advanced LyoSmart operating system, assisting users to optimize the freeze-drying process.

- Display of ice condenser temperature and vacuum degree
- Product temperature display according to sublimation pressure curve
- Manual or automatic mode process optional
- Programmable freeze-drying steps with high flexibility
- Vacuum degree of each step can be set and controlled independently
- Historical data recorder for convenient analysis of processing
- USB port to download data
- Password access protection



► **Full vacuum degree control**

- The vacuum degree and duration of each step can be set and controlled at auto mode
- Vacuum control to reduce drying times by up to 40%
- Real-time vacuum degree is displayed on the screen
- Pirani vacuum gauge with accuracy 0.001 mbar



► **Various “one button” mode**

- One-button start: in auto mode, press "START" to start process according to preset program
- One-button stop: vacuum valve shut off, preventing product polluted due to vacuum back suction of pump oil
- One-button vent: aeration valve controls vacuum vent gradually, preventing product polluted due to vacuum back suction of pump oil
- One-button defrosting: hot air defrosting of cold trap, to enter a new process faster

► **Advanced cooling system**

- Built-in ice condenser coils for higher efficiency
- -60°C: single compressor system, high efficiency
- -90°C: 2-compressor cascade system, cold trap temperature reaches -90°C within 15 minutes
- SECOP brand compressor for durable use

► **Optimized flange seal**

- Flat surface & L-shaped sealing ring design, no need for extra manual pressure
- Vacuum grease not required, easier to use and clean
- Proven excellent sealing performance for long-term use



► **Choose the system fits your application best!**

Product designation format:

LyoPro **4/8/12** **S/M**
 Product series Ice capacity (in kg) Temperature level

► **Ice capacity**

Max. ice capacity	Performance	Models
4 kg	4 kg/24h	LyoPro 4S, LyoPro 4M
8 kg	6 kg/24h	LyoPro 8S, LyoPro 8M
12 kg	8 kg/24h	LyoPro 12S, LyoPro 12M

► **Temperature level**

Temperature level	Temperature	Cooling system	Typical application area
S	-60°C	1 compressor	Aqueous products
M	-90°C	2 compressors cascade	Products with low freezing points

*For organic solvent application, please choose LyOrganic series freeze dryers.

► **The UNIVERSAL package with special price**

Start your lyophilization with our universal package, which offers both manifolds and shelves, fulfilling most applications for successful processes every day.

► **The LyoPro series UNIVERSAL package contains**

- LyoPro freeze dryer main machine
- Pirani vacuum gauge
- Pressure control valve
- Acrylic chamber with 8 rubber connection valves
- Rack supporting board
- 250ml & 500ml flask (qty refer to specification sheet)
- 4-layer shelves with unique ShelFlex™ design
- Vacuum pump
- Vacuum hose



► **ShelFlex™ design**

Unique ShelFlex™ design, various configurations with different number of shelves and spacing are available.

Shelves				
Configuration	Number	Φ	Area	Spacing
1	2	265 mm	0.10 m ²	140 mm
2	3	265 mm	0.15 m ²	105 mm
3	4	265 mm	0.20 m ²	70 mm
4	8	265 mm	0.40 m ²	35 mm

*For configuration 4, another 4 pcs of shelves to be ordered as optional.

► **Advantages**

- All-inclusive package, avoid confusion in accessory selection
- ShelFlex™ design, adjustable shelf number and spacing with flexibility
- Cost-effective, fulfill more applications with lower cost
- Fast delivery, UNIVERSAL package always in stock for urgent demand

► **Various freeze-drying packages for different applications**



Manifold		Shelves				
No.	Number of connection valves	Number	Φ	Area	Spacing	Stopper
①	8	4	265 mm	0.20 m ²	70 mm	-
②	8	2	265 mm	0.10 m ²	140 mm	●
③	-	4	265 mm	0.20 m ²	70 mm	-
④	-	2	265 mm	0.10 m ²	140 mm	●
⑤	8	-	-	-	-	-
⑥	12	-	-	-	-	-
⑦	24	-	-	-	-	-

*Package 1, also named as "universal package", is recommended and supplied as standard configuration.

● Basic equipment - Not available



Specification (-60 °C)

Model		LyoPro 4S	LyoPro 8S	LyoPro 12S
Basic	Type	Desktop	Desktop	Floor-standing
Ice condenser	Temperature	-60°C	-60°C	-60°C
	Max. ice capacity	4 kg	8 kg	12 kg
	Chamber volume	6.5 L	15 L	21 L
	Performance	4 kg/24h	6 kg/24h	8 kg/24h
Control	Display	7" touch screen	7" touch screen	7" touch screen
	Operating system	LyoSmart	LyoSmart	LyoSmart
	Ice condenser temperature display	●	●	●
	Product temperature sublimation	●	●	●
	Vacuum degree display	●	●	●
	Vacuum gauge type	Pirani	Pirani	Pirani
	Ultimate vacuum (mbar)	0.020	0.020	0.020
	Vacuum venting	Automatic	Automatic	Automatic
	Controlled gas (N ₂ , clean air, etc) backfill	●	●	●
	Hot gas defrosting	●	●	●
Data	Data record	●	●	●
	USB port	●	●	●
	RS232 port	○	○	○
Cooling system	Compressor	SECOP compressor x 1	SECOP compressor x 1	SECOP compressor x 1
Power	Voltage (V/Hz)	220/50;220/60;120/60	220/50;220/60;120/60	220/50;220/60;120/60
	Rated power	800 W	900 W	1000 W
Noise level	Freeze dryer	48 (dB)	52 (dB)	52 (dB)
Dimensions	Main unit (W x D x H)	455*550*430 mm	800*533*400 mm	700*650*1015 mm
Weight	Main unit	53 kg	69 kg	125 kg
Universal drying package	Drying chamber	Acrylic	Acrylic	Acrylic
	Rack	4-shelf	4-shelf	4-shelf
	Rack supporting board	●	●	●
	Number of rubber connection valves	8	8	8
	Round bottom flask*	250ml*6 & 500ml*6	250ml*6 & 500ml*6	250ml * 6 & 500ml * 6
Vacuum pump	Vacuum pump	●	●	●
	Oil mist filter	●	●	●
	Vacuum pump hose	●	●	●

*: Flasks are supplied free of charge. Compensation is not applicable in case of damage during transportation.

The data provided refers to the base unit with ambient conditions of +10 °C to +25 °C. Subject to change without prior notice. ● Basic equipment ○ Optional

 **Specification (-90 °C)**


Model		LyoPro 4M	LyoPro 8M	LyoPro 12M	
Basic	Type	Desktop	Desktop	Floor-standing	
Ice condenser	Temperature	-90°C	-90°C	-90°C	
	Max. ice capacity	4 kg	8 kg	12 kg	
	Chamber volume	6.5 L	15 L	21 L	
	Performance	4 kg/24h	6 kg/24h	8 kg/24h	
Control	Display	7" touch screen	7" touch screen	7" touch screen	
	Operating system	LyoSmart	LyoSmart	LyoSmart	
	Ice condenser temperature display	●	●	●	
	Product temperature sublimation	●	●	●	
	Vacuum degree display	●	●	●	
	Vacuum gauge type	Pirani	Pirani	Pirani	
	Ultimate vacuum (mbar)	0.015	0.015	0.015	
	Vacuum venting	Automatic	Automatic	Automatic	
	Controlled gas (N2, clean air, etc) backfill	●	●	●	
	Hot gas defrosting	●	●	●	
	Data	Data record	●	●	●
		USB port	●	●	●
RS232 port		○	○	○	
Cooling system	Compressor	SECOP compressor x 2	SECOP compressor x 2	SECOP compressor x 2	
Power	Voltage (V/Hz)	220/50	220/50;220/60;120/60	220/50;220/60;120/60	
	Rated power	1300 W	1600 W	1800 W	
Noise level	Freeze dryer	48 (dB)	52 (dB)	52 (dB)	
Dimensions	Main unit (W x D x H)	455*550*430 mm	800*533*400 mm	700*650*1015 mm	
Weight	Main unit	66 kg	82 kg	125 kg	
Universal drying package	Drying chamber	Acrylic	Acrylic	Acrylic	
	Rack	4-shelf	4-shelf	4-shelf	
	Rack supporting board	●	●	●	
	Number of rubber connection valves	8	8	8	
	Round bottom flask*	250ml*6 & 500ml*6	250ml*6 & 500ml*6	250ml * 6 & 500ml * 6	
Vacuum pump	Vacuum pump	●	●	●	
	Oil mist filter	●	●	●	
	Vacuum pump hose	●	●	●	

*: Flasks are supplied free of charge. Compensation is not applicable in case of damage during transportation.

The data provided refers to the base unit with ambient conditions of +10 °C to +25 °C. Subject to change without prior notice. ● Basic equipment ○ Optional

LyOrganic series

Freeze dryer for organic solvent

LyOrganic series is designed for applications dealing with both organic solvents and inorganic solvents. Typical applications include organic dyes and pigments, polymers in benzol, organic substances in solvents, preparative HPLC fractions, solvent removal after chromatography, etc.

Condenser chamber temperature reaches -105°C , vacuum degree is controlled during the process, oil-free vacuum pump is supplied as standard, meeting high requirement of organic solvent drying.

LyOrganic series freeze dryers are compact in design, high-efficiency in performance, smart in operation, and dependable in daily use.



► Resistant to organic solvents

Typical examples include HPLC fractions with organic or inorganic solvents, such as acetonitrile, TFA and other alcohols, or other products with t-butanol, DMSO, etc.

► Special system configurations for safety

LyOrganic series freeze dryers are compact in design, high-efficiency in performance, smart in operation, and dependable in daily use.

- Durable: high-quality materials for excellent chemical resistance, such as PTFE coated stainless steel ice condenser, and solvent-resistant seals
- High performance internal condenser, most solvents are frozen or liquefied in the chamber at -105°C

► Innovative features

- Compact, high-performance laboratory systems with a small footprint
- 7" touch screen installed with LyoSmart operating system
- Programmable freeze-drying steps
- Controlled vacuum degree for optimized process times
- Built-in ice condenser coils for higher efficiency
- Fast defrost via integrated hot gas
- Automatic non-return venting valve
- Low noise level: 54 dB(A)

▶ LyoSmart operating system

7" touch screen controller, installed with advanced LyoSmart operating system, assisting users to optimize the freeze-drying process.

- Display of ice condenser temperature and vacuum
- Indirect product temperature determination based on vapor pressure curve
- Manual or automatic drying process optional
- Programmable freeze-drying steps
- Vacuum degree of each step can be set and controlled independently
- Historical data recorder for convenient analysis of processing
- USB port to download data
- Password access protection



▶ Vacuum degree control

- The vacuum degree and duration of each step can be set and controlled at auto mode
- Vacuum control to reduce drying times by up to 40%
- Real-time vacuum degree is displayed on the screen

▶ Various "one button" mode

- One-button start: in auto mode, press "START" to start process according to preset program
- One-button stop: vacuum valve shut off automatically
- One-button vent: aeration valve controls vacuum vent gradually
- One-button defrosting: hot air defrosting of cold trap, to enter a new process faster

▶ Advanced cooling system

- Built-in ice condenser coils for higher efficiency
- -105°C: 2-compressor cascade system, suitable for organic solvents with ultra low freezing points
- SECOP brand compressor for durable use

▶ Ease of use and durability

- Flat surface & L-shaped sealing ring design, no need for extra manual pressure
- Vacuum grease not required, easier to use and clean
- Proven excellent sealing performance for long-term use
- PTFE coated stainless steel ice condenser, resistant to organic solvents



▶ Choose the system fits your application best!

Product designation format

LyOrganic **8/12** **L** **Drying Package**
 Product series Ice capacity (in kg) Temperature level T/U

▶ The various models have different maximum ice capacities:

Max. ice capacity	Performance	Ice condenser temperature	Models
8 kg	6 kg/24h	-105°C	LyOrganic 8LT, LyOrganic 8LU
12 kg	8 kg/24h	-105°C	LyOrganic 12LT, LyOrganic 12LU

The LyOrganic series are available in two defined packages, fully aligned to the requirements of freeze drying with organic solvents.

► Organic solvent package – manifold

- LyOrganic freeze dryer main machine
- Drying manifold with 12 chemical-resistant connection valves for flasks
- 250ml & 500ml flask (qty refer to specification sheet)
- Vacuum sensor
- Vacuum control to reduce drying time by up to 40%
- Drain valve
- Internal condenser, stainless steel with solvent-resistant PTFE material coating
- High performance oil-free vacuum pump, resistant to organic solvents
- Vacuum pump hose



► Organic solvent package – manifold and shelves

- LyOrganic freeze dryer main machine
- Drying chamber with 8 rubber connection valves
- Rack supporting board
- 250ml & 500ml flask (qty refer to specification sheet)
- 4-layer shelves with unique ShelfFlex™ design
- Vacuum sensor
- Vacuum control to reduce drying time by up to 40%
- Drain valve
- Internal condenser, stainless steel with solvent-resistant PTFE material coating
- Solvent-resistant seals and connectors
- High performance oil-free vacuum pump, resistant to organic solvents
- Vacuum pump hose



► ShelfFlex™ design

Unique ShelfFlex™ design, various configurations with different number of shelves and spacing are available.

Shelves				
Configuration	Number	Φ	Area	Spacing
1	2	265 mm	0.10 m ²	140 mm
2	3	265 mm	0.15 m ²	105 mm
3	4	265 mm	0.20 m ²	70 mm
4	8	265 mm	0.40 m ²	35 mm

*For configuration 4, another 4 pcs of shelves to be ordered as optional.

► Advantages

- All-inclusive package, avoid confusion in accessory selection
- ShelfFlex™ design, adjustable shelf number and spacing with flexibility
- Cost-effective, fulfill more applications with lower cost
- Fast delivery, UNIVERSAL package always in stock for urgent demand

► Various freeze-drying packages for different applications



Manifold		Shelves				
No.	Number of connection valves	Number	Φ	Area	Spacing	Stopper
①	8	4	265 mm	0.20 m ²	70 mm	-
②	8	2	265 mm	0.10 m ²	140 mm	●
③	-	4	265 mm	0.20 m ²	70 mm	-
④	-	2	265 mm	0.10 m ²	140 mm	●
⑤	8	-	-	-	-	-
⑥	12	-	-	-	-	-
⑦	24	-	-	-	-	-

● Basic equipment - Not available

 **Specification**


Model		LyOrganic 8LT	LyOrganic 8LU	LyOrganic 12LT	LyOrganic 12LU
Basic	Type	Desktop	Desktop	Floor-standing	Floor-standing
Ice condenser	Temperature	-105°C	-105°C	-105°C	-105°C
	Max. ice capacity	8 kg	8 kg	12 kg	12 kg
	Chamber volume	15 L	15 L	21 L	21 L
	Performance	6 kg/24h	6 kg/24h	8 kg/24h	8 kg/24h
Control	Display	7" touch screen	7" touch screen	7" touch screen	7" touch screen
	Operating system	LyoSmart	LyoSmart	LyoSmart	LyoSmart
	Ice condenser temperature display	●	●	●	●
	Product temperature sublimation	●	●	●	●
	Vacuum degree display	●	●	●	●
	Vacuum gauge	●	●	●	●
	Ultimate vacuum (mbar)*	0.015	0.015	0.015	0.015
	Vacuum venting	Automatic	Automatic	Automatic	Automatic
	Controlled gas (N ₂ , clean air, etc) backfill	●	●	●	●
	Hot gas defrosting	●	●	●	●
Data	Data record	●	●	●	●
	USB port	●	●	●	●
	RS232 port	○	○	○	○
Cooling system	Compressor	SECOP compressor x 2			
Power	Voltage (V/Hz)	220/50;220/60;120/60	220/50;220/60;120/60	220/50;220/60;120/60	220/50;220/60;120/60
	Rated power	1700 W	1700 W	1900 W	1900 W
Noise level	Freeze dryer	52 (dB)	52 (dB)	52 (dB)	52 (dB)
Dimensions	Main unit (W x D x H)	800*533*400 mm	800*533*400 mm	700*650*1015 mm	700*650*1015 mm
Weight	Main unit	83 kg	83 kg	145 kg	145 kg
Manifold package for flasks	Manifold type	T-type	-	T-type	-
	Number of vessel ports	12	-	12	-
	Round bottom flask**	250ml * 6 & 500ml *6	-	250ml * 6 & 500ml *6	-
Universal drying package	Drying chamber	-	Acrylic	-	Acrylic
	Rack	-	4-shelf	-	4-shelf
	Rack supporting board	-	●	-	●
	Number of vessel ports	-	8	-	8
	Round bottom flask*	-	250ml * 6 & 500ml *6	-	250ml * 6 & 500ml * 6
Vacuum pump	Vacuum pump SVF-100DR (oil-free)	●	●	●	●
	Vacuum pump NDD12 (with oil)	○	○	○	○
	Vacuum pump hose	●	●	●	●

*: The ultimate vacuum 0.015 mbar is achieved when working with vacuum pump NDD12 (NAVAC brand).

● Basic equipment ○ Optional - Not available

** : Flasks are supplied free of charge. Compensation is not applicable in case of damage during transportation.

The data provided refers to the base unit with ambient conditions of +10 °C to +25 °C. Subject to change without prior notice.

Pilot Freeze Dryer, 2-chamber system

LyoPilot D Series

Equipped with stainless steel shelf areas with synthetic silicone heat transfer fluid, the larger 8/12kg pilot systems meet the highest standards of the pharma and biotech industries. The 2-chamber system is integrated with advanced PAT functions for optimized process.

LyoPilot D series includes temperature control system, vacuum level control system, identical operation and freeze-drying process controls, are ideally suited for freeze-drying solids or liquids in a wide variety of containers, such as vials, flasks, trays, microtiter plates, and ampoules. LyoPilot D series is extremely valuable to optimize (shorten) the length of the freeze-drying cycle, especially if there is potential for process repetition or scale-up for production.



Advanced system technology


► 2-chamber systems

- Optimal vapor transport.
- Standard side-opening ice condenser for easier observation and cleaning.
- Ultra-low pressure drop between the two chambers.
- Excellent drying rates for sensitive materials and products with low eutectic points or low collapse temperatures.
- Intermediate valve for determining the transition from primary to secondary drying phase.

► Choose the system fits your application best!

Product designation format:

LyoPilot **D** **8/12** **S/-**
 Product series Dual (2) chamber Ice capacity (kg) With/without stoppering

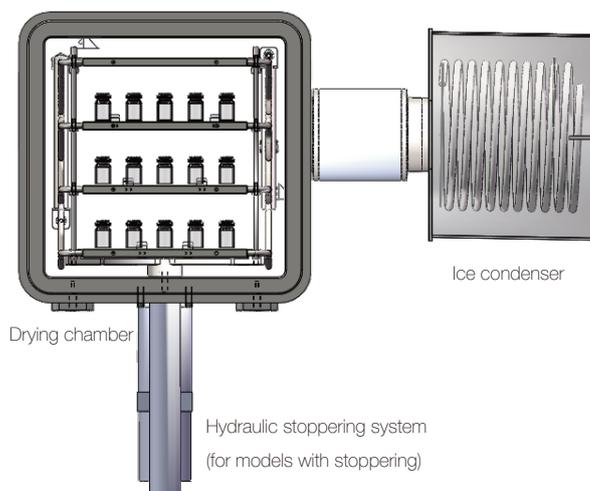

► Silicone heat transfer fluid shelves

- Equipped with stainless steel shelf areas with silicone heat transfer fluid, providing excellent thermal conductivity.
- The accuracy of different layers and a same layer is within 1°C.
- Hydraulic stopper system, rising from bottom to top and easy for cleaning (available for LyoPilot D8S/LyoPilot D12S)

Design features

▶ Minimal space requirements

Compact design with highly efficient and economical operation.

▶ Precise temperature and vacuum control

- Precise Temperature and Vacuum Control Thanks to independently developed adaptive PID algorithms, the equipment achieves higher precision in temperature and vacuum control.
- The 10mm thick stainless steel silicone oil shelf ensures a temperature difference of 1°C.
- The ASV adaptive gas injection system, achieving a vacuum control accuracy of up to 0.001 mbar.

▶ Comprehensive PAT (Process Analytical Technology) tools

Equipped with advanced tools, providing clearer and more intuitive control over the freeze-drying process:

- Eutectic point testing
- Resistivity monitoring
- Pressure rise test
- Pressure comparison method
- Flash nucleation technology

▶ Outstanding refrigeration system

- Ice condenser: 2-compressor cascade system, ice condenser temperature reaches -90°C.
- Shelf: pulse-type electronic expansion valve, featuring a wide regulation range, quick and precise control on shelf temperature, and low energy consumption.

▶ Efficient & convenient ice condenser

- Internal ice condenser, condenser coil contacts vapor directly, for easier vapor capturing and lower energy consumption.
- Standard hot gas defrosting, enables customers to quickly proceed to the next freeze-drying process.
- Mechanism door lock of ice condenser, for convenient cleaning.

▶ Stopping function (available for LyoPilot D8S/LyoPilot D12S)


The fully hydraulic automatic stopping system operates from bottom to top, saving time and facilitating cleaning.

Unique system controller LyoSMART+

LyoSMART+ system controller make an easy-to-use and intuitive user interface. All of the extensive accessories are also integrated. Reproducible results are assured by automatic process sequences.

► Convenient and intuitive

- Color touchscreen with clear display.
- Automatic or manual mode optional.
- Intuitive program entry, using various freeze-drying sequences and recipes.
- Memory space for 30 user-defined programs with 36 steps.
- Real-time freezing drying diagram display, selective parameter display, data query & export.
- The main interface is clear and concise, facilitating program editing and operation.
- The control system can set safety values such as vacuum level and resistivity for process feedback control.
- Password protection for startup and lock screen protection when operator is away during processing.
- Equipped with USB port & RS232 port.

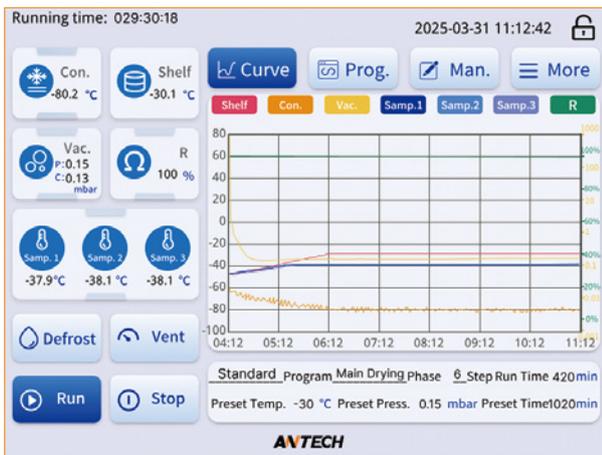


Phase	Step	TEMP. (°C)	Vac. (mbar)	Heating/Cooling Time Settings (min)	Duration (min)
Pre-freezing	1	4	--	10	60
Pre-freezing	2	-40	--	20	60
Pre-freezing	3	-40	--	10	180
Pump preheat	4	--	--	--	--
Main Drying	5	-35	0.15	10	30
Main Drying	6	-20	0.15	60	600
Main Drying	7	0	0.15	60	600
Final Drying	8	25	0.05	60	360

Smart solutions for optimal process observation

► Freezing point determination

- The LyoTR sensor temperature & resistivity combined monitors both electrical resistance and product temperature. From the curves of both of these variables, you can determine the eutectic point and freezing point of your product.
- Our LyoSOFT software and LyoSMART+ makes this possible. The LyoTR sensor allows automated control of the energy supply to the shelves during the main drying phase, so you can avoid critical temperatures during the main drying phase. This reduces the risk of defrosting effects on the product.

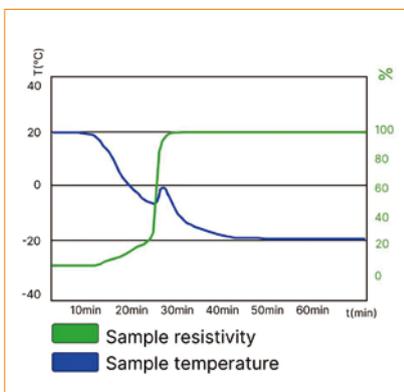


Vial supporter (optional)

- The ultra-fine probes with flexible wire harness is compact, highly resistant to interference, corrosion-resistant, and easy to calibrate.

► Pressure rise test

- The pressure rise test is achieved through the design of the intermediate valve between the drying chamber and the ice condenser. During the drying process, this valve is closed to measure the pressure increase in the drying chamber over a defined time period, which serves as a key indicator for determining the freezing point.
- Typically, the pressure rise rate at the end of main drying should be within 0.01 mbar/min. The accurate vacuum control is critical for this function. Thanks to Antech's ultra-high accurate ASV (Advanced Smart Valve) vacuum control, our pressure rise test delivers superior reliability.

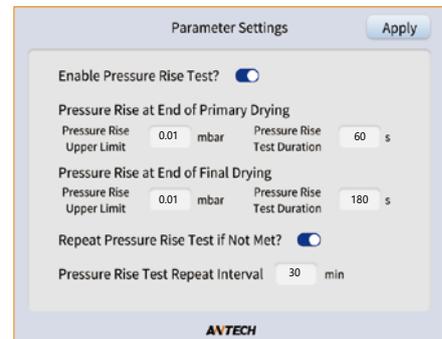


► Pre-vacuum and micro-negative pressure stoppering

- The range accuracy of the traditional Pirani vacuum gauge is limited. Pirani data is not available at higher ranges (10-1000mbar). We can add a silicon crystal vacuum gauge to supplement the measurement and control of this part of the range.
- After putting the sample into the machine, you can first set a rough vacuum degree of 600-950mbar for the airtightness check of the equipment to ensure that the machine is in a good airtight condition and to prevent airtightness problems caused by too low temperature, which leads to vacuum abnormalities.
- After the freeze-drying is completed, a higher vacuum degree of 500-800mbar can also be set. The equipment can backfill the inert gas according to the set value, so that the vial can be plugged under a slightly negative pressure, thereby ensuring that the dried sample can be stored for a longer time.

► Comparative pressure measurement

- The end of the main drying phase can be detected by using two different vacuum measurement sensors (Pirani gauge and capacitive sensor).
- When the difference between the pressure measurements falls below a preselected threshold, final drying is started automatically.

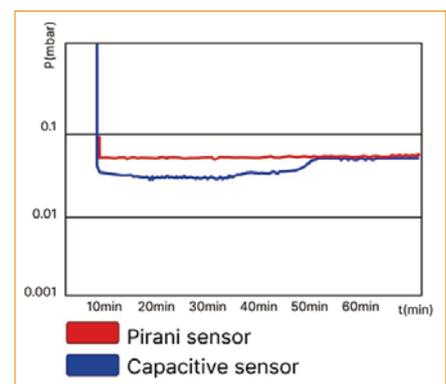

► Instantaneous Co-Crystallization system

- Due to thermodynamic factors, the crystallization process of samples within the same batch is not synchronized but random. Variations in supercooling time lead to differences in crystal structure, resulting in inconsistent freeze-drying outcomes (e.g., moisture residue, appearance, etc.).
- The instantaneous crystallization technology can integrate vacuum changes in the cold trap chamber and sample chamber, enabling sudden condensation of samples through controlled vacuum adjustments. Unlike other ice fog-induced nucleation methods, our instantaneous crystallization technology does not require additional auxiliary devices or exogenous substances, even when handling small-scale sample tests.

LyoSOFT: operate on PC



Time	Shelf	Com	Vac1	Vac2	Samp.1	Samp.2	Samp.3	R	Program	
1	2025-03-20 11:25:51	17.4	18.8	1000	1.287	16.7	16.8	16.7	100	FreeD1_1
2	2025-03-20 11:25:56	17.4	18.8	1000	1.287	16.7	16.7	16.8	100	FreeD1_1
3	2025-03-20 11:26:01	17.3	18.8	1000	1.287	16.7	16.8	16.8	100	FreeD1_1
4	2025-03-20 11:26:06	17.4	18.8	1000	1.287	16.7	16.7	16.7	100	FreeD1_1
5	2025-03-20 11:28:11	17.4	18.8	1000	1.287	16.7	16.8	16.7	100	FreeD1_1
6	2025-03-20 11:28:16	17.5	18.8	1000	1.287	16.7	16.8	16.7	100	FreeD1_1
7	2025-03-20 11:28:21	17.5	18.7	1000	1.287	16.7	16.8	16.7	100	FreeD1_1
8	2025-03-20 11:28:26	17.6	18.7	1000	1.287	16.7	16.7	16.7	100	FreeD1_1
9	2025-03-20 11:28:31	17.6	18.8	1000	1.287	16.7	16.8	16.8	100	FreeD1_1
10	2025-03-20 11:28:36	17.6	18.7	1000	1.287	16.7	16.8	16.7	100	FreeD1_1
11	2025-03-20 11:28:41	17.7	18.7	1000	1.287	16.7	16.8	16.7	100	FreeD1_1
12	2025-03-20 11:28:46	17.7	18.8	1000	1.287	16.7	16.7	16.7	100	FreeD1_1
13	2025-03-20 11:28:51	17.7	18.7	1000	1.287	16.7	16.8	16.7	100	FreeD1_1
14	2025-03-20 11:28:56	17.6	18.7	1000	1.287	16.7	16.8	16.7	100	FreeD1_1
15	2025-03-20 11:29:01	17.6	18.7	1000	1.287	16.7	16.8	16.7	100	FreeD1_1
16	2025-03-20 11:29:06	17.6	18.7	1000	1.287	16.7	16.7	16.7	100	FreeD1_1
17	2025-03-20 11:29:11	17.6	18.7	1000	1.287	16.7	16.8	16.7	100	FreeD1_1
18	2025-03-20 11:29:16	17.9	18.7	1000	1.287	16.7	16.7	16.7	100	FreeD1_1
19	2025-03-20 11:29:21	17.9	18.7	1000	1.287	16.7	16.8	16.7	100	FreeD1_1
20	2025-03-20 11:29:26	17.8	18.7	1000	1.287	16.7	16.8	16.8	100	FreeD1_1

► For precise documentation and evaluation

- Our experience shows that freeze-drying processes must be precisely monitored and documented. Regardless of the drying parameters or batch size, this is the only way to ensure accurate analysis.
- The optional LyoSOFT software (PC based) enables comprehensive recording and archiving of all process data. Installed on a separate computer, the software receives data directly from the freeze dryer via an RS232 interface.
- With LyoSOFT, programs for freeze drying can be developed and process data can be viewed in real time in graphical format. Operation is consistent and uniform across all unit sizes, as LyoSOFT is also used with larger production freeze drying systems.

► Are you planning to scale up?

- Data recording on USB drive.
- Simple process documentation with LyoSOFT.
- LyoSOFT for process control and documentation.
- Process monitoring with the LyoTR sensor to avoid undesired defrosting effects.
- Automatic determination of the freezing point for reliable process control.
- High-precision temperature and vacuum control ensures the accuracy of the freeze-drying process.
- Scientific structural layout not only meets the convenience of operation, but also ensures the repeatability of the process.
- LyoSOFT software complies with current GAMP standards.

With stoppering

► Specification

Model		LyoPilot D8S	LyoPilot D12S
Basic	Type	2-chamber	2-chamber
	Stopper	●	●
Ice condenser	Type	Internal placed	Internal placed
	Temperature *	-90°C	-90°C
	Max. ice capacity	8 kg	12 kg
	Chamber volume	18L	30L
	Performance	6kg/24h	10kg/24h
	Hot gas Defrosting	●	●
Shelf	Type	Silicone heat transfer fluid	Silicone heat transfer fluid
	Dimensions(W*D)	240*410mm	280*500mm
	Number of shelves	3+1	4+1
	Area	0.3m ²	0.56m ²
	Spacing	80mm	73mm
	Hydraulic stoppering system	Bottom to top	Bottom to top
	Temperature range	-55°C ~ +60°C	-60°C ~ +60°C
	Temperature accuracy	±1°C	±1°C
	Temperature uniformity	±1°C	±1°C
	Cooling rate (+20°C to -40°C)	1.3°C/min	1.3°C/min
Controller	Display	8" touch screen	8" touch screen
	Operating system	LyoSmart+	LyoSmart+
	Password protection	●	●
Data	Data record	●	●
	USB port	●	●
	RS232 port	●	●
Vacuum	Vacuum gauge type	Pirani type x 1, capacitive type x 1	Pirani type x 1, capacitive type x 1
	Ultimate vacuum	0.008mbar	0.008mbar
	Vacuum control accuracy	0.001mbar	0.001mbar
	Vacuum leakage rate	0.005mbar L / s	0.005mbar L / s
Process control	Pressure rise test	●	●
	Safety pressure	●	●
	T&R monitoring to prevent defrosting	●	●
	Target/actual comparison T _{product} and T _{shelf} area	●	●
PAT tools	T&R(determining freezing point, product resistance)	●	●
	Product temperature measurement	●	●
	Capacitive pressure measurement	●	●
	Comparative pressure measurement	●	●
Backfill Aeration	Automatic backfill (N ₂ , clean air, etc)	●	●
	Automatic aeration	●	●
Cooling	Cooling system	2 Compressor cascade	2 Compressor cascade
	Refrigerant	Hydrocarbon	Hydrocarbon
Power	Voltage	220V/50Hz;220V/60Hz;120V/60Hz	380V/50Hz
	Rated power	3400W	4500W
Noise level	Freeze dryer	56 (dB)	60 (dB)
Dimensions	Freeze dryer (W x H x D)	960*765*1320mm	1065*860*1380mm
Weight	Freeze dryer	330Kg	450Kg
Vacuum pump	Vacuum pump and hose	●	●
Software for PC	LyoSOFT	○	○

*: Tested at ambient temperature 24°C, with good ventilation

● Basic configuration ○ Optional

Without stoppering

► Specification

Model		LyoPilot D8	LyoPilot D12
Basic	Type	2-chamber	2-chamber
	Stopper	/	/
Ice condenser	Type	Internal placed	Internal placed
	Temperature *	-90°C	-90°C
	Max. ice capacity	8 kg	12 kg
	Chamber volume	18L	30L
	Performance	6kg/24h	10kg/24h
	Hot gas Defrosting	●	●
Shelf	Type	Silicone heat transfer fluid	Silicone heat transfer fluid
	Dimensions(W*D)	315*410mm	350*500mm
	Number of shelves	3+1	4+1
	Area	0.39m ²	0.7m ²
	Spacing	80mm	73mm
	Temperature range	-55°C ~ +60°C	-60°C ~ +60°C
	Temperature accuracy	±1°C	±1°C
	Temperature uniformity	±1°C	±1°C
	Cooling rate (+20°C to -40°C)	1.3°C/min	1.3°C/min
Controller	Display	8" touch screen	8" touch screen
	Operating system	LyoSmart+	LyoSmart+
	Password protection	●	●
Data	Data record	●	●
	USB port	●	●
	RS232 port	●	●
Vacuum	Vacuum gauge type	Pirani type x 1, capacitive type x 1	Pirani type x 1, capacitive type x 1
	Ultimate vacuum	0.008mbar	0.008mbar
	Vacuum control accuracy	0.001mbar	0.001mbar
	Vacuum leakage rate	0.005mbar L / s	0.005mbar L / s
Process control	Pressure rise test	●	●
	Safety pressure	●	●
	T&R monitoring to prevent defrosting	●	●
	Target/actual comparison T _{product} and T _{shelf} area	●	●
PAT tools	T&R(determining freezing point, product resistance)	●	●
	Product temperature measurement	●	●
	Capacitive pressure measurement	●	●
	Comparative pressure measurement	●	●
Backfill Aeration	Automatic backfill (N ₂ , clean air, etc)	●	●
	Automatic aeration	●	●
Cooling	Cooling system	2 Compressor cascade	2 Compressor cascade
	Refrigerant	Hydrocarbon	Hydrocarbon
Power	Voltage	220V/50Hz;220V/60Hz;120V/60Hz	380V/50Hz
	Rated power	3400W	4500W
Noise level	Freeze dryer	56 (dB)	60 (dB)
Dimensions	Freeze dryer (W x H x D)	960*765*1320mm	1065*860*1380mm
Weight	Freeze dryer	300Kg	420Kg
Vacuum pump	Vacuum pump and hose	●	●
Software for PC	LyoSOFT	○	○

*: Tested at ambient temperature 24°C, with good ventilation

● Basic configuration ○ Optional

► Specification

LyoPilot D8S (with stoppering)

Shelf dimensions (W x D): 240x 410 mm								
Vial volume (total)			2ml	6ml	10ml	20ml	50ml	100ml
Number of shelves	Area (m ²)	Spacing (mm)	Max. number of vials *					
3(Standard)	0.30	80	1263	630	540	360	N/A	N/A
2(Customized)*	0.20	120	842	420	360	240	100	64

LyoPilot D12S (with stoppering)

Shelf dimensions (W x D): 280x 500 mm								
Vial volume (total)			2ml	6ml	10ml	20ml	50ml	100ml
Number of shelves	Area (m ²)	Spacing (mm)	Max. number of vials *					
4(Standard)	0.56	73	2380	1200	1008	648	N/A	N/A
2(Adjustable)*	0.28	146	1190	600	504	324	160	100

LyoPilot D8 (without stoppering)

Shelf dimensions (W x D): 315x 410 mm								
Vial volume (total)			2ml	6ml	10ml	20ml	50ml	100ml
Number of shelves	Area (m ²)	Spacing (mm)	Max. number of vials *					
3(Standard)	0.39	80	1611	852	714	429	N/A	N/A
2(Customized)*	0.26	120	1074	568	476	286	144	100

LyoPilot D12 (without stoppering)

Shelf dimensions (W x D): 350x 500 mm								
Vial volume (total)			2ml	6ml	10ml	20ml	50ml	100ml
Number of shelves	Area (m ²)	Spacing (mm)	Max. number of vials *					
4(Standard)	0.70	73	2952	1508	1296	800	N/A	N/A
2(Adjustable)*	0.35	146	1476	754	648	400	196	122

*Adjustable: equipment is adjustable at user's site based on standard configuration.

*Customized: non-standard configuration to be manufactured at factory.

*Note: Data for maximum load, qty will be less when using loading frames.

Vial size overview

Total vial volume	2 ml	6 ml	10 ml	20 ml	50 ml	100 ml
Vial type	2R	6R	10R	20R	50H	100H
↓ mm vial only	35	40	45	55	73	95
↓ mm with Lyo plug	45	50	55	65	83	105
φ mm	16	22	24	30	43	52
Net fill volume at 1 cm fill height (ml)	1.2	2.2	4.0	4.6	6.0	7.0



Qingdao Antech Scientific Co., Ltd.

Tel: +86 532 87890321

Email: info@antechscientific.com

Web: www.antechscientific.com



TABLE OF CONTENTS

▶ **LyoPro series** Lab freeze dryer 1-5

▶ **LyOrganic series** Lab freeze dryer for organic solvent 6-9

▶ **LyoPilot D Series** Pilot freeze dryer, 2-chamber system 10-17

