

# LabSafe Series


## -25°/MDF-25U270/400/600

LabSafe series upright lab freezers are used to store medicines, vaccines or lab materials in pharmacies, hospitals, clinics and labs. Product temperature sensor is put into aluminum block to simulate real product temperature. Access port is available for additional temperature probe.

These units use environmentally friendly, natural refrigerants that are fully compliant with new EPA/SNAP/UL regulations.

Temperature data storage and USB port are supplied as standard. A pdf format file including temperature curve and temperature data can be downloaded. A quick check of maximum and minimum temperature record in previous 24 hours is available, making the job much easier for operators.

▼ MDF-25U270

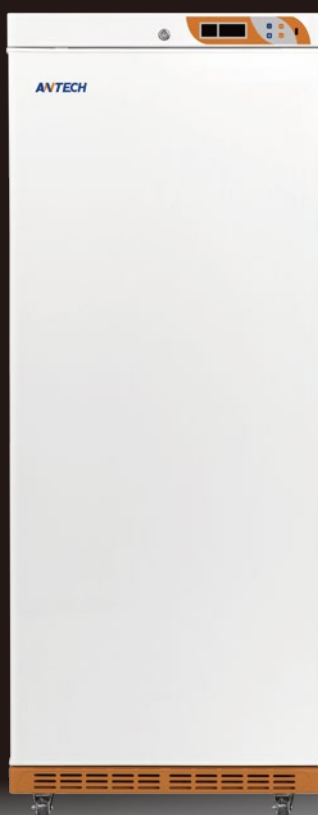

▼ MDF-25U600



NATURAL  
REFRIGERANT



USB PORT



HIGH/LOW  
TEMPERATURE ALARM



POWER FAILURE  
ALARM



DOOR AJAR A  
LARM



REVERSIBLE  
DOORS

MODEL	MDF-25U270	MDF-25U400	MDF-25U600
Temperature range (°C)	-10 ~ -25	-10 ~ -25	-10 ~ -25
Capacity (L)	280	380	550
External dimensions W×D×H (mm)	600×560×1540	600×650×1805	775×750×1905
Shelves (pcs)	5	6	6
Data storage & USB port	Yes	Yes	Yes
Power supply	220V/50Hz 220V/60Hz 120V/60Hz	220V/50Hz 220V/60Hz 120V/60Hz	220V/50Hz 220V/60Hz 120V/60Hz