

PharmaStore series pharmacy refrigerators

Antech PharmaStore series pharmacy refrigerators are designed to store drugs, vaccine and reagents in drug stores, clinics, hospitals and laboratories.

Advanced forced air cooling & dual sensor control systems are adopted for excellent temperature uniformity and accurate temperature control.

PharmaStore series are equipped with complete professional functions to guarantee safety of products, such as high/low temperature alarm, door ajar alarm, sensor error alarm, power failure alarm, remote alarm contact, etc.

Available with a wide range of capacities, from 60 liter desktop refrigerator to 1015 liter side-by-side door refrigerator. PharmaStore series refrigerators are compliant with main international requirements, including CE, UL, etc.

This units uses environmentally friendly, natural refrigerant fully compliant with new EPA/SNAP/UL regulations.

TC-110 Controller

Microprocessor electronic controller for PharmaStore series pharmacy refrigerators, advanced dual sensor control system for accurate temperature control.



LOW CO₂
EMISSIONS



NATURAL
REFRIGERANT



CFC /HCFC
FREE



HIGH/LOW
TEMPERATURE ALARM



POWER FAILURE
ALARM



DOOR AJAR A
LARM



SAFETY DOOR
LOCK



FORCE AIR
COOLING



CE
CERTIFICATION

PERFORMANCE

- Forced air cooling for excellent temperature uniformity
- Variable (cycle) defrost
- Microprocessor controller with LED digital display
- Dual sensor design: Independent control sensor and display sensor

ERGONOMIC DESIGN

- Interior LED lighting for convenient operation
- Ultra-quiet design reduces noise level
- Adjustable feet for leveling

SAFETY

- High/low temperature alarm, door ajar alarm, sensor error alarm, power failure alarm with backup battery
- Safety door lock, preventing unauthorized access
- Access port for additional temperature probe
- Remote alarm contact

ENVIRONMENT FRIENDLY

- Natural refrigerant compliant with EPA and SNAP non-ozone depleting mandates
- High efficiency cooling system for lower energy consumption and carbon emission



MPR-60



MPR-120



MPR-280



MPR-362




MPR-440




- ▶ **Safety door lock**
Preventing unauthorized access



- ▶ **Forced air cooling**
Excellent temperature uniformity and fast recovery after door open

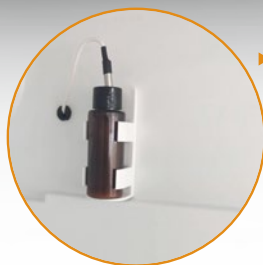

- ▶ **Dual sensor system**
Product sensor in Emery sand bottle to simulate real product temperature; Control sensor in air for accurate temperature control.

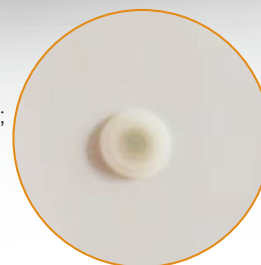

- ▶ **Access port**
For easy access of additional temperature probe



- ▶ **Remote alarm contact**



- ▶ **Backup battery**
For alarm in case of power failure

Specification

Model		MPR-60	MPR-120	MPR-280	MPR-362	MPR-440
Main Specifications	Temperature Range (℃)	2～8	2～8	2～8	2～8	2～8
	Climate Class	N	N	N	N	N
	Capacity (L)	60	120	280	362	440
	Cooling Type	Forced-Air	Forced-Air	Forced-Air	Forced-Air	Forced-Air
	Defrost	Auto	Auto	Auto	Auto	Auto
	Refrigerant	R600a	R600a	R600a	R600a	R600a
	Noise (dB)	45	45	45	45	45
Control System	Controller	TC-110	TC-110	TC-110	TC-110	TC-110
	Display	LED digital	LED digital	LED digital	LED digital	LED digital
	Display Accuracy	0.1 ℃	0.1 ℃	0.1 ℃	0.1 ℃	0.1 ℃
	Sensor	Display sensor & control sensor				
Power	Power Supply (V/Hz)	220/50, 220/60, 120/60				
	Power (W)	65	80	135	240	280
Dimensions	Inner (W*D*H, mm)	394*376*418	477*424*632	473*550*1160	473*550*1460	473*550*1760
	External (W*D*H, mm)	480*535*580	597*540*880	595*710*1340	595*710*1640	595*710*1940
	Packing (W*D*H, mm)	520*565*620	635*585*935	640*760*1370	640*760*1670	640*760*1970
	Net Weight (KG)	22	37	53	59	65
	Gross Weight (KG)	26	45	64	71	78
Alarm Functions	Visible/audible Alarm	Yes	Yes	Yes	Yes	Yes
	High/Low Temperature	Yes	Yes	Yes	Yes	Yes
	Sensor Error	Yes	Yes	Yes	Yes	Yes
	Door Ajar	Yes	Yes	Yes	Yes	Yes
	Power Failure	Yes	Yes	Yes	Yes	Yes
	Low Battery	Yes	Yes	Yes	Yes	Yes
	Remote Alarm Contact	Yes	Yes	Yes	Yes	Yes
Materials	Exterior	Painted steel	Painted steel	Painted steel	Painted steel	Painted steel
	Interior	Polystyrene	Polystyrene	Polystyrene	Polystyrene	Polystyrene
	Insulation	Polyurethane	Polyurethane	Polyurethane	Polyurethane	Polyurethane
	Door	Glass door x 1	Glass door x 1	Glass door x 1	Glass door x 1	Glass door x 1
	Inner light	LED	LED	LED	LED	LED
	Shelf	Coated Steel Wire	Coated Steel Wire	Coated Steel Wire	Coated Steel Wire	Coated Steel Wire
	Shelf Qty	1 large & 1 small	3 large & 1 small	5 large & 1 small	7 large & 1 small	9 large & 1 small
Others	Door Lock	Yes	Yes	Yes	Yes	Yes
	Test Access Port	Yes	Yes	Yes	Yes	Yes
	Adjustable Feet	Yes	Yes	Yes	Yes	Yes
	Backup battery	Yes	Yes	Yes	Yes	Yes