

LyoPro series

Laboratory freeze dryer

Optimize and simplify your freeze-drying process

Antech Scientific offers advanced LyoPro series freeze-drying system with unique LyoSmart operating system, realizing full control of vacuum degree, making the freeze-drying process simpler, faster and better.

LyoPro series freeze dryers are compact in design, high-efficiency in performance, smart in operation, and in daily use.



► Innovative features

- Compact, high-performance laboratory systems with a small footprint
- 7" touch screen installed with LyoSmart operating system
- Programmable freeze-drying steps
- Controlled vacuum degree for optimized process times
- Built-in ice condenser coils for higher efficiency
- Fast defrost via integrated hot gas
- Automatic non-return venting valve, prevents "vacuum pump" oil mist from entering cold trap
- Pirani vacuum sensor with accuracy 0.001 mbar
- Low noise level: 48~52dB(A)
- High quality vacuum pump

► LyoSmart operating system


7" touch screen controller, installed with advanced LyoSmart operating system, assisting users to optimize the freeze-drying process.

- Display of ice condenser temperature and vacuum degree
- Product temperature display according to sublimation pressure curve
- Manual or automatic mode process optional
- Programmable freeze-drying steps with high flexibility
- Vacuum degree of each step can be set and controlled independently
- Historical data recorder for convenient analysis of processing
- Password access protection



► Full vacuum degree control

- The vacuum degree and duration of each step can be set and controlled at auto mode
- Vacuum control to reduce drying times by up to 40%
- Real-time vacuum degree is displayed on the screen
- Pirani vacuum gauge with accuracy 0.001 mbar

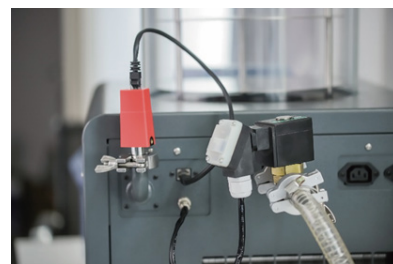

► Various “one button” mode

- One-button start: in auto mode, press "START" to start process according to preset program
- One-button stop: vacuum valve shut off, preventing product polluted due to vacuum back suction of pump oil
- One-button vent: aeration valve controls vacuum vent gradually, preventing product polluted due to vacuum back suction of pump oil
- One-button defrosting: hot air defrosting of cold trap, to enter a new process faster

► Advanced cooling system

- -55°C: single compressor system, high efficiency
- -85°C: 2-compressor cascade system, cold trap temperature reaches -85°C within 10 minutes
- SECOP brand compressor for durable use

► Optimized flange seal

- Flat surface & L-shaped sealing ring design, no need for extra manual pressure
- Vacuum grease not required, easier to use and clean
- Proven excellent sealing performance for long-term use



► Choose the system fits your application best!

Product designation format:

LyoPro **4/8/12** **S/M**
 Product series Ice capacity (in kg) Temperature level

► Ice capacity

Max. ice capacity	Performance	Models
4 kg	4 kg/24h	LyoPro 4S, LyoPro 4M
8 kg	6 kg/24h	LyoPro 8S, LyoPro 8M
12 kg	8 kg/24h	LyoPro 12S, LyoPro 12M

► Temperature level

Temperature level	Temperature	Cooling system	Typical application area
S	-55°C	1 compressor	Aqueous products
M	-85°C	2 compressors cascade	Products with low freezing points

*For organic solvent application, please choose LyOrganic series freeze dryers.

► The UNIVERSAL package with special price

Start your lyophilization with our universal package, which offers both manifolds and shelves, fulfilling most applications for successful processes every day.

► The LyoPro series UNIVERSAL package contains:

- LyoPro freeze dryer main machine
- Pirani vacuum gauge
- Pressure control valve
- Acrylic chamber with 8 rubber connection valves
- 250ml & 500ml flask (qty refer to specification sheet)
- 4-layer shelves with unique ShelfFlex™ design
- Vacuum pump
- Vacuum hose



► ShelfFlex™ design

Unique ShelfFlex™ design, various configurations with different number of shelves and spacing are available.

Shelves				
Configuration	Number	Φ	Area	Spacing
1	2	265 mm	0.10 m ²	140 mm
2	3	265 mm	0.15 m ²	105 mm
3	4	265 mm	0.20 m ²	70 mm
4	8	265 mm	0.40 m ²	35 mm

*For configuration 4, another 4 pcs of shelves to be ordered as optional.

► Advantages

- All-inclusive package, avoid confusion in accessory selection
- ShelfFlex™ design, adjustable shelf number and spacing with flexibility
- Cost-effective, fulfill more applications with lower cost
- Fast delivery, UNIVERSAL package always in stock for urgent demand

► Various freeze-drying packages for different applications



Manifold		Shelves				
No.	Number of connection valves	Number	Φ	Area	Spacing	Stopper
①	8	4	265 mm	0.20 m ²	70 mm	-
②	8	2	265 mm	0.10 m ²	140 mm	●
③	-	4	265 mm	0.20 m ²	70 mm	-
④	-	2	265 mm	0.10 m ²	140 mm	●
⑤	8	-	-	-	-	-
⑥	12	-	-	-	-	-
⑦	24	-	-	-	-	-

*Package 1, also named as "universal package", is recommended and supplied as standard configuration.

● Basic equipment

- Not available



Specification (-55°C)

Model		LyoPro 4S	LyoPro 8S	LyoPro 12S
Basic	Type	Desktop	Desktop	Floor-standing
Ice condenser	Temperature	-55°C	-55°C	-55°C
	Max. ice capacity	4 kg	8 kg	12 kg
	Chamber volume	6.5 L	15 L	21 L
	Performance	4 kg/24h	6 kg/24h	8 kg/24h
Control	Display	7" touch screen	7" touch screen	7" touch screen
	Operating system	LyoSmart	LyoSmart	LyoSmart
	Ice condenser temperature display	●	●	●
	Product temperature sublimation	●	●	●
	Vacuum degree display	●	●	●
	Vacuum gauge type	Pirani	Pirani	Pirani
	Ultimate vacuum (mbar)	0.02	0.02	0.02
	Vacuum control range (mbar)	0.001~1000	0.001~1000	0.001~1000
	Vacuum venting	Automatic	Automatic	Automatic
	Controlled gas (N ₂ , clean air, etc) backfill	●	●	●
	Hot gas defrosting	●	●	●
Data	Data record	●	●	●
	USB port	○	○	○
	RS232 port	○	○	○
Cooling system	Number of compressor	1	1	1
	Compressor brand	SECOP	SECOP	SECOP
Power	Voltage (V/Hz)	220/50;220/60;120/60	220/50;220/60;120/60	220/50;220/60;120/60
	Rated power	800 W	900 W	1000 W
Noise level	Freeze dryer	48 (dB)	52 (dB)	52 (dB)
Dimensions	Main unit (W x D x H)	455*550*430 mm	800*533*400 mm	700*650*1015 mm
Weight	Main unit	53 kg	69 kg	125 kg
Universal drying package	Drying chamber	Acrylic	Acrylic	Acrylic
	Shelves	4	4	4
	Number of rubber connection valves	8	8	8
	Round bottom flask*	250ml*6 & 500ml*6	250ml*6 & 500ml*6	250ml * 6 & 500ml * 6
Vacuum pump	Vacuum pump	●	●	●
	Oil mist filter	●	●	●
	Vacuum pump hose	●	●	●

*: Flasks are supplied free of charge. Compensation is not applicable in case of damage during transportation.

The data provided refers to the base unit with ambient conditions of +10 °C to +25 °C. Subject to change without prior notice. ● Basic equipment ○ Optional

► **Specification (-85 °C)**



Model		LyoPro 4M	LyoPro 8M	LyoPro 12M
Basic	Type	Desktop	Desktop	Floor-standing
Ice condenser	Temperature	-85°C	-85°C	-85°C
	Max. ice capacity	4 kg	8 kg	12 kg
	Chamber volume	6.5 L	15 L	21 L
	Performance	4 kg/24h	6 kg/24h	8 kg/24h
Control	Display	7" touch screen	7" touch screen	7" touch screen
	Operating system	LyoSmart	LyoSmart	LyoSmart
	Ice condenser temperature display	●	●	●
	Product temperature sublimation	●	●	●
	Vacuum degree display	●	●	●
	Vacuum gauge type	Pirani	Pirani	Pirani
	Ultimate vacuum (mbar)	0.01	0.01	0.01
	Vacuum control range (mbar)	0.001~1000	0.001~1000	0.001~1000
	Vacuum venting	Automatic	Automatic	Automatic
	Controlled gas (N2, clean air, etc) backfill	●	●	●
	Hot gas defrosting	●	●	●
Data	Data record	●	●	●
	USB port	○	○	○
	RS232 port	○	○	○
Cooling system	Number of compressor	2	2	2
	Compressor brand	SECOP	SECOP	SECOP
Power	Voltage (V/Hz)	220/50	220/50;220/60;120/60	220/50;220/60;120/60
	Rated power	1300 W	1600 W	1800 W
Noise level	Freeze dryer	48 (dB)	52 (dB)	52 (dB)
Dimensions	Main unit (W x D x H)	455*550*430 mm	800*533*400 mm	700*650*1015 mm
Weight	Main unit	66 kg	82 kg	125 kg
Universal drying package	Drying chamber	Acrylic	Acrylic	Acrylic
	Shelves	4	4	4
	Number of rubber connection valves	8	8	8
	Round bottom flask*	250ml*6 & 500ml*6	250ml*6 & 500ml*6	250ml * 6 & 500ml * 6
Vacuum pump	Vacuum pump	●	●	●
	Oil mist filter	●	●	●
	Vacuum pump hose	●	●	●

*: Flasks are supplied free of charge. Compensation is not applicable in case of damage during transportation.

The data provided refers to the base unit with ambient conditions of +10 °C to +25 °C. Subject to change without prior notice. ● Basic equipment ○ Optional